

GB 2000 Series

2 Channels GAS Mixer



Each Channel: 0 – 5000 mL/min

Accuracy: 1.0%

Repeatability: 0,10% of reading

Response time: 500 ms

Software bundle.

For STANDARD Flows

Total Mixture Flow Rate: Up to 10 L/min



Up to 2 Channels.

GB 2000 can manage and mix up to 2 different sources of gas with incredible high precision, creating accurate gas mixtures.

Touch Screen.

Our Touch Display allows you to manage and control any value and parameter of your flows in real time and without effort.

Compatible.

The integrability with ProfiBus, ModBus, & Matlab Simulink gives the possibility of an easy integration into other systems.

The Ultimate Gas Mixers.


MCQ Instruments offers the best, cheapest and fastest solutions to manage, control, measure and create Gas Mixtures or Dilutions.

We produce advanced, certified, interconnected, and trusted technologies that can help you to achieve your goals without any efforts. From laboratories to industrial needs, we have the right product that will manage your gas processes with automation and an incomparable precision.

Gases:

N₂,
O₂,
CO₂,
CH₄, Air,
He, H₂

(additional gases supported on request)

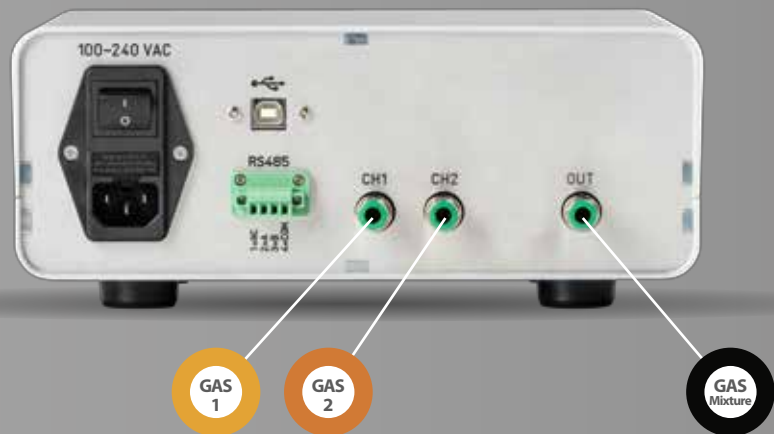

Touch Display.

The Gas Blender 2000 Series is a stand-alone device, and it allows you to control all the parameters and access all the main features by the Touch Display. The user interface already contains several of the most important actions such as Play and Stop the gas mixing, regulation of the gas ratio for each channel, native gases setting and much more.

Easy Connected.

Gas Blender 2000 Series is not only USB connected, but Bluetooth ready, with the possibility of future Wi-Fi implementation.

This functionality allows customers to remote control the instrument. Its integrability with ProfiBus, ModBus, and Matlab Simulink gives the possibility of an extremely easy integration into other systems.



Intuitive Software.



MCQ Mixture Creator is our mixing software bundled for free, that allows direct access and control to all of our products. Thanks to a user-friendly interface, you can manage the gas mixing process by intervening directly on the single channel to modify and control both the final blending and your dynamic dilution.

Certified Products.



Pending

Pending



Automated Feedback System.

With Gas Blender 2000 Series, you can create, save and load your automatic mixture program, set the percentage of each channel, the total flow of the gas mixture and processing duration.

The program code can be edited anytime, exported to an external file and reloaded later.

Market Analysis.

Research and Development (R&D) Labs, Institutes and Centers

- Calibration of Instruments
- Preparation of Test Gases
- Analytical Equipment
- Local Preparation of Gas Mixtures
- Chemical Laboratory

Analysis

- Gas Chromatography
- Quality Control

Food & Beverage

- Regulation of Controlled Atmospheres

Life Sciences


- Biotechnology
- Regulation of Gaseous Atmospheres
- Studies on Cell Cultures
- Pharmaceutical

Accurate. Compact. Reliable.


“The MCQ Gas Blender 2000 Series is the ultimate solution specifically designed for laboratories, for non-aggressive gas mixings and dilutions. It is an innovative all-in-one replacement for traditional mass flow controllers multi-channel solutions.”

Specifications.


Accuracy: (for each channel)	10-100% FS: 0,2% of FS < 10% FS: 1% of FS
Repeatability:	0,10% of reading
Response time: (for each channel)	500 ms
Power Supply:	In Bundle
Operation Pressure:	Max 3 bar
Working Temperature:	0-50 °C
Mass Flow Rates: (for each channel)	0 – 5000 mL/min (standard)
Total Mixture Flow Rate:	Up to 10 L/min
Gases:	N2, O2, CO2, CH4, Air, He, H2 (additional gases supported on request)
Compatibility:	Profibus, Modbus, Matlab Simulink
Communication:	USB PC Interface RS485 with open proprietary protocol
Fittings:	Swagelok NPT or Push-in fittings for 6-4mm tubes. Other on request
Inlet Gas Channel:	2 Inputs
Outlet Gas Channel:	1 Output of mixed gases
User Interface:	Gas mixing software. Touch Display.
Certifications:	CE, TUV, cTUVUS, Made in Italy

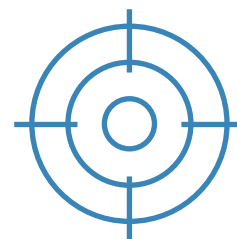


High Accuracy
Native calibration on available gas with 1% of accuracy for each channel

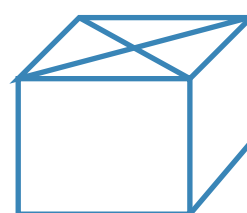

Lab in a Box
All the tools in a single, compact and easy to transport box.

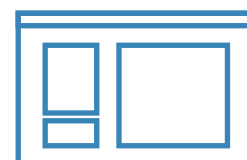

Software Bundle
Mixture Creator is an user friendly interface, with which you can easily manage any gas mixing process.

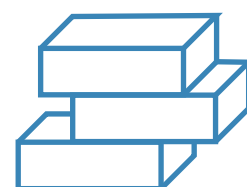

GB Family
We have as many solutions as many are your needs.